PermaQuick®

Bonding System

Description:

1. PermaQuick Primer is a tissue kind, ethyl alcohol, one part primer which contains natural resins that aid penetration of hydrophilic substrates such as dentin. Cured PermaQuick Primer is somewhat elastic. It provides the optimum preparation for the PermaQuick Bonding Resin.
2. PermaQuick Bonding Resin is a 45% filled, fluoride releasing, light cured bonding resin that penetrates into primer treated tubules and inter-tubular spaces. PermaQuick Bonding Resin is as radiopaque as dentin.

Indications:

Use PermaQuick for most bonding needs in restorative dentistry. PermaQuick is conducive for bonding to Dentin/enamel/composite, using composite or amalgam and Porcelain Metal Light cure luting.

Procedure:

Remove luer lock caps and attach the appropriate tip: UltraEtch: Inspiral Brush tip (some prefer the Blue Micro tip). PermaQuick Primer: Black Mini Brush tip. PermaQuick Bonding Resin: Inspiral Brush tip.

DENTIN/ENAMEL/COMPOSITE

1. Surface Preparation: Prepared dentin/enamel/composite surfaces to be bonded should be clean and isolated from contamination of soft tissues and oral fluids. A rubber dam (DermaDam®) is recommended if possible.
2. Etch: Apply Ultra-Etch® onto preparation surfaces and leave for 15 seconds. Thoroughly wash to remove etch. Blow or vacuum off excess water, leaving surface slightly wet.
3. Apply a puddled coat of PermaQuick Primer. Air thin primer to evaporate the alcohol solvent carrier (approximately 12 seconds). Light cure 20 seconds using standard halogen light or 10 seconds using high energy light.
4. Apply a puddled coat of PermaQuick bonding resin onto primed surface. Gently air thin. Do not overthin (no dry spots, maintain a glossy surface). Light cure 20 seconds using standard halogen light or 10 seconds using high energy light.
5. Restore:
 - a. For Composite: We recommend Vit-I-escence™ or Amelogen® composites. Build with incremental layering keeping the first layer thin (0.5 mm). Followed by 2-3 mm increments. Sculpt and shape. Polish to a luster and seal with PermaSeal® composite sealer.
 - b. For Amalgam: Follow steps 1-4 and condense mixed amalgam into cavity preparation. Sculpt and shape appropriately.

Note: Highest bond strengths occur when the surface has been left visibly moist; however, respectable bond strength occurs even when the dentin has been significantly dried.

PORCELAIN

Remove luer lock caps and attach the appropriate tips:

Ultradent® Porcelain Etch: Inspirial Brush tip. Ultradent® Silane: Black Mini Brush tip

Important Note: Before clinical procedure, dispense Porcelain Etch onto paper pad to verify and ensure evenness of flow. If there is resistance to flow, replace tip with a new one and verify flow before proceeding.

1. Surface Preparation: Remove loose and/or thin projections of porcelain and bevel as needed. Freshen fractured surface of porcelain with medium diamond bur and/or micro-abrade.
2. Isolate: Isolate adjacent tissues and restorations before etching. When rubber dam placement is not possible, isolation can be established with OpalDam® and EtchArrest®.
3. Etch: Apply Ultradent Porcelain Etch to prepared surface. Etch 1 minute and neutralize acid with EtchArrest until bubbling stops. Rinse and dry. Etched surface should be frosty. If not, repeat Porcelain Etch step.
4. Silanate: Brush Ultradent Silane onto etched porcelain. Leave for 1 minute, followed by air dry to evaporate any remaining carrier.
5. Apply puddled coat of PermaQuick bonding resin onto silanated surface. Gently air thin. Do not over-thin. Light cure 20 seconds using a standard halogen light or 10 seconds using a high energy light.
6. Restore: Use appropriate composite, building incrementally. We recommend Vit-I-escence or Amelogen. Build with incremental layering. Keeping the first layer thin (0.5mm) followed by 2-3mm increments. Note: Highest bond values are achieved with (1) Preparation of porcelain surface, (2) Ultradent Porcelain Etch, (3) Ultradent Silane, and (4) PermaQuick Bonding Resin (without using PermaQuick Primer). PermaQuick Primer application to silanated surface provides high, respectable bond values but does not provide highest test values obtained.

METAL

1. Surface Preparation: Micro-abrade non-precious metal surface.
Note: Bond strength to metal is determined by metal alloy and surface preparation technique as follows (highest to lowest bond values):
 - a. Precious metal: Micro-abrade, then tin plate.
 - b. Non-precious metal: Micro-abrasion is best.
 - c. Medium diamond bur preparation is better than smooth finish, but not as good as A or B
2. Clean: In all cases (except if tin plating), metal surface should be etched with Ultra-Etch for 15 seconds, followed by rinsing thoroughly for 5 seconds with firm air/water spray.
3. Prime: (Priming is necessary only if dentin is involved.) Apply a puddled coat of primer onto the dentin surface. Following application, use air to evaporate the alcohol solvent carrier. Light cure for 20 seconds.
4. Bonding Resin: Apply a puddled coat of PermaQuick Bonding Resin into primed surface for 15 seconds. Gently air thin. Do not over-thin (no dry spots or bare dentin). Light cure 20 seconds.
5. Restore: Restorative of choice may now be placed. We recommend Vit-I-escence or Amelogen. Start with 0.5mm initial layer. Cure and build layers in 2-3mm increments.

LUTING

1. Surface Preparation Tooth Structure: Isolate preparation with displacement cord and appropriate hemostatic agent. Clean preparation thoroughly.
2. Etch: Etch/condition the preparation 15 seconds with Ultra-Etch followed by a 5 second air/water wash. Remove excess water, leaving surface slightly wet.
3. Prime: Apply a puddled coat of PermaQuick Primer into preparation. Following application, gently air dry 5- 10 seconds to evaporate the alcohol solvent carrier. Light cure 20 seconds.
4. Bonding Resin: Apply a puddled coat of PermaQuick Bonding Resin onto primed surface. Gently air thin. Do not over thin. Do not light cure.
Note: Complete polymerization of both the bonding resin and the light cured luting resin occurs after placing veneer and initiating light cure. This ensures complete seating of veneer.
5. Restore: Express 0.5mm layer of PermaFlo® flowable composite (or when chemical cure is necessary PermaFlo DC) into appropriately etched and silanated porcelain veneer. Using the all plastic Black Mini tip will prevent leaving a dark metal abraded line on the rough porcelain. Seat veneer, removing excess resin, and light tack to hold in place. Apply DeOx® oxygen barrier solution to prevent formation of oxygen inhibition layer. Light cure 40-60 seconds from the lingual. Clean excess and polish margins. Do not use more than one thin layer of die spacer. This will minimize dimensional change upon polymerization shrinkage. If more than one thin layer of die spacer is used, the luting resin will pull the porcelain in toward the tooth upon polymerization rather than separating from the dentin, as can occur with weaker bonding agents

Precautions:

1. Remember to establish adequate isolation (rubber dam, etc.) before applying Ultra-Etch or Porcelain Etch. Isolation is also important to prevent contamination during primer and bonding resin application. Use EtchArrest to protect tissues and neutralize hydrofluoric acid for situations in which a rubber dam is not feasible.
2. It is imperative to experiment with etchant application on something inert such as a paper pad BEFORE using on patients. This will allow you to become familiar with the ease of delivery and the viscous properties of the material. If flow does not come easily, stop and replace tip.
3. Confine phosphoric and hydrofluoric acids to the areas being etched or cleaned. Remember that these acids can irritate and burn tissues. Be careful not to splash into eyes, onto mucous membranes, skin or any other areas. Any tissues or surfaces inadvertently exposed to acid should be immediately rinsed with lots of water.
4. Any resin can be sensitizing. Wash exposed area completely to remove the material. If dermatitis, rash or any other allergic type reaction is observed, consult a physician.
5. Remember that light activated bonding resins are sensitive to ambient light; syringe caps should be replaced following use to prevent light exposure.
6. Resins should be kept in refrigerator if not used on a daily basis. Room temperature is adequate for storage of product used daily.
7. Replace luer lock caps when syringes are not in use. TipSoc covers may be used on tips for storage of syringes from day to day.
8. When dispensing intraorally, never force syringe plungers. If there is excess resistance to flow, replace tip with a new one and verify flow before proceeding. Use only recommended tips with each specific material.

PermaQuick

Bonding System

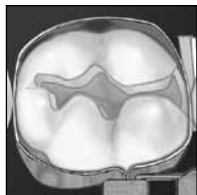
Figure 1: Dentin/Enamel or Composite: Rubber dam placed, cavity prep completed.

Figure 2: Apply Ultra-Etch® 15 seconds. Air/water wash 5 seconds. Remove surface water, do not desiccate.

Figure 3: Apply a puddled coat of PermaQuick Primer. Air thin 12 seconds. Light Cure 20 seconds with standard light or 10 seconds with high energy light.

Figure 4: Apply a puddled coat of PermaQuick Bonding onto primed surface. Gently air thin. Light cure 20 seconds with standard light or 10 seconds with high energy light.

Porcelain Repair Kit Precautions:

1. Remember to establish adequate isolation (rubber dam, etc.) before applying Ultra-Etch or Porcelain Etch. Isolation is also important to prevent contamination during primer and bonding resin application. Use EtchArrest to protect tissues and neutralize hydrofluoric acid for situations in which a rubber dam is not feasible.
2. It is imperative to experiment with etchant application on something inert such as a paper pad BEFORE using on patients. This will allow you to become familiar with the ease of delivery and the viscous properties of the material. If flow does not come easily, stop and replace tip.
3. Confine phosphoric and hydrofluoric acids to the areas being etched or cleaned. Remember that these acids can irritate and burn tissues. Be careful not to splash into eyes, onto mucous membranes, skin or any other areas. Any tissues or surfaces inadvertently exposed to acid should be immediately rinsed with lots of water.
4. Any resin can be sensitizing. Wash exposed area completely to remove the material. If dermatitis, rash or any other allergic type reaction is observed, consult a physician.
5. Remember that light activated bonding resins are sensitive to ambient light; syringe caps should be replaced following use to prevent light exposure.
6. Resins should be kept in refrigerator if not used on a daily basis. Room temperature is adequate for storage of product used daily.

Figure 5: Restore tooth using restorative material of choice.

Figure 6: Porcelain: Apply Ultradent® Porcelain etch 1 minute. Rinse and dry. Apply Ultradent® Silane; air dry 1 minute. Bond as show in Figures 1 thru 5.

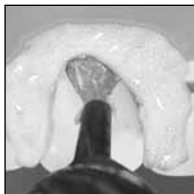
Figure 7: Cast metal: Microabrade (tin plate precious metals). Rinse and dry. Bond as show in Figures 1 thru 5.

Figure 8: Luting: Follow etch, prime and bond technique shown in Figures 1 thru 4 with exception of light curing. Cure bond resin and luting resin simultaneously.

7. Replace luer lock caps when syringes are not in use. TipSoc covers may be used on tips for storage of syringes from day to day.
8. When dispensing intraorally, never force syringe plungers. If there is excess resistance to flow, replace tip with a new one and verify flow before proceeding. Use only recommended tips with each specific material.

General Precautions:

1. For Professional use only.
2. Review instructions, precautions, and MSDS before beginning treatment. Use only as directed.
3. Keep products out of heat/sunlight.
4. Avoid skin exposure to resins.
5. Isolate strong chemicals to area of treatment.
6. Confirm that patient has no known allergies to treatment materials.
7. Test flow of materials from syringe and tip before using intraorally.
8. Never force syringe plungers.
9. Use Ultradent® Syringe Covers and/or clean and disinfect syringes between patients.

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

PermaQuick Primer

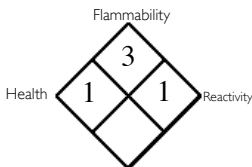For product MSDS please see our website: www.ultradent.com

PermaQuick Bonding Resin

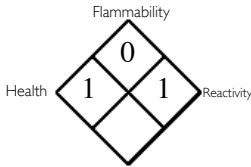For product MSDS please see our website: www.ultradent.com

Hazard Rating
4 = Severe
3 = Serious
2 = Moderate
1 = Slight
0 = Minimal

Ultradent syringes have an expiration date stamped on the side of the syringe consisting of one letter and three numbers. The letter is a lot number used for manufacturing purposes and the three numbers are the expiration date. The first two numbers are the month, and the third number is the last number of the year.

U.S. federal law restricts the sale of this device by or on the order of a state licensed dental professional.
Keep out of reach of children.

For immediate reorder and/or complete descriptions of Ultradent's product line, refer to Ultradent's catalog or call Toll Free 1-800-552-5512.

Outside U.S. call (801) 572-4200
or visit www.ultradent.com.

© Copyright 2005 Ultradent Products, Inc. U.S. Patent Nos.: 4,997,371; 5,269,684; 5,534,562; 5,700,148; 5,708,052 and 5,722,833.. Other U.S. and International Patents Granted and Pending. All Rights Reserved.